



**AQUATEC® OCEAN  
ERGO SHOWER & COMMODE CHAIRS**



**Yes, you can.®**

**AQUATEC®**

# Customer insight-led design to enhance independence, safety and flexibility.

Personal care is a highly private activity and is an important part of a person's life and wellbeing. The **Aquatec Ocean Ergo Family** of shower and commode chairs have been designed with user comfort, independence and dignity in mind. Based on the well-established Aquatec Ocean line, the Ocean Ergo family comes packed with a range of new features, designed to make personal care activities easier for both the user and caregiver.

## **Safety, independence and dignity**

A newly contoured seat surface has been developed to encourage a more upright and stable seating position, promoting user independence and overall

comfort. Enhanced stability at the front of the base allows for a smaller footprint and offers users added safety and reassurance when seated in the chair. Furthermore, the optimized wheelbase makes it even easier for a caregiver to maneuver the chair in small spaces as well as allowing for a better fit over standard toilets for independent toileting.

## **Modular and flexible design**

The modular design of the Ocean Ergo line means that the chairs can be adapted to suit a wide range of user needs by making some simple adjustments, or by adding some of our easy to fit, off-the-shelf accessories.

## **OTHER features include:**

- ▶ Stable, stainless steel frame
- ▶ Machine washable, tension adjustable back with quick-release buckles
- ▶ Easy disassembly for transportation and storage



# Features

## 1 Innovative ergonomic seat

- ▶ Insight-led design encourages a more upright seating position for enhanced independence and comfort

## 2 Greater adaptability

- ▶ With an increased height adjustment range on the footrests and armrests

## 4 Compact and easy to handle

- ▶ Optimized wheel base makes it easy to maneuver in small spaces

## 3 Extremely stable

- ▶ Enhanced stability at the front for additional safety and reassurance





## Contoured, ergonomic seat plate

- ▶ Based on the imprints people leave behind when seated
- ▶ A 5° seat tilt lifts the knees into a natural squat position for toileting
- ▶ Transfer supports at the front to allow the user to position themselves into a comfortable position



## Highly Configurable

Available with an extensive range of positioning, safety, toileting and seating accessories that can be used to support a wide range of conditions and postural needs.


See full range of accessories online at [www.invacare.com](http://www.invacare.com)



**Positioning**



**Safety**

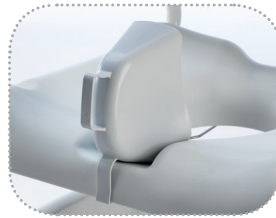

**Toileting**



**Seating and Backrests**

## Tilt-in-space models:

**Ocean Ergo VIP and Ocean Ergo Dual VIP**

- ▶ Both have a -5° to 40° seat tilt angle with a smart weight shift mechanism that ensures optimal stability on a small footprint, while the Ocean Dual VIP has an additional backrest recline of 5° to 35°



# Introducing the Ocean Ergo Family

From a standard transit model, to tilt and recline, meet the flexible line of shower chairs that offer personal care solutions for a wide range of physical conditions.

## Modular personal care solutions

SAFE | INDEPENDENT | ERGONOMIC



### **Ocean 24” Ergo XL ▲**

An entry level, plus-size model with a maximum user weight of 400 lbs

### **Ocean Ergo XL ▲**

An entry level plus-size attendant propelled model, with a maximum user weight of 400 lbs

### **Ocean VIP Ergo ▲**

Designed for individuals who require slightly more assistance with postural control when seated

### **Ocean Ergo ▲**

An entry level attendant propelled model

### **Ocean Dual VIP Ergo ▲**

Designed for individuals with complex seating and positioning needs

### **Ocean 24” Ergo ▲**

An entry level self-propelling model

## Technical data



Seat width ▼



Seat depth ▼



Total width ▼



Total depth ▼

Width between armrests ▼

<b>Ocean Ergo</b>	19"/480 mm	19"/480 mm	22"/565 mm	35.5"/905 mm	18"/455 mm
<b>Ocean Ergo XL</b>	19"/480 mm	19"/480 mm	25"/645 mm	35.5"/905 mm	21"/535 mm
<b>Ocean 24" Ergo</b>	19"/480 mm	19"/480 mm	26.5"/670 mm	42"/1050 mm	18"/455 mm
<b>Ocean 24" Ergo XL</b>	19"/480 mm	19"/480 mm	26.5"/670 mm	42"/1050 mm	21"/535 mm
<b>Ocean VIP Ergo</b>	19"/480 mm	19"/480 mm	22"/565 mm	39.5"/1000 mm	18"/455 mm
<b>Ocean Dual VIP Ergo</b>	19"/480 mm	19"/480 mm	22"/565 mm	39.5"/1000 mm	18"/455 mm



Tilt angle seat ▼



Tilt angle back ▼




Weight ▼

Load capacity ▼

<b>Ocean Ergo</b>	-	-	37.5 lbs/17 kg	330 lbs/150 kg
<b>Ocean Ergo XL</b>	-	-	40 lbs/18 kg	400 lbs/180 kg
<b>Ocean 24" Ergo</b>	-	-	51 lbs/23 kg	330 lbs/150 kg
<b>Ocean 24" Ergo XL</b>	-	-	53 lbs/24 kg	400 lbs/180 kg
<b>Ocean VIP Ergo</b>	-5° to 40°	-	52 lbs/23.5 kg	330 lbs/150 kg
<b>Ocean Dual VIP Ergo</b>	-5° to 40°	5° to 35°	62.5 lbs/28.3 kg	330 lbs/150 kg

## Frame Color

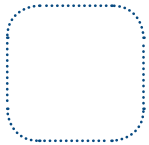

White

## Regulatory Data



\*All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further details.

\*\*Excluding tires and wheels

NOTE: Rollators are not intended to be used as wheelchairs

**Invacare Corporation**  
www.invacare.com

**USA**  
One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

**Canada**  
570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2019  
Invacare Corporation.  
All rights reserved.  
Trademarks are identified  
by the symbols ™ and  
®. All trademarks are  
owned by or licensed  
to Invacare Corporation  
unless otherwise noted or  
its subsidiaries. Aquatec  
is owned by Invacare Int'l  
GmbH. Specifications are  
subject to change without  
notification.  
Form No. 19-116  
190807



**Yes, you can.®**