



**Yes, you can.®**

## Sling Specifications

	Divided-leg				Full-body			Full-body w/Commode			Toileting
	P	M	L	XL	M	L	XL	M	L	XL	L
<b>Dimensions</b>											
Width	33"	37.5"	41"	44.5"	41.5"	45.5"	45.5"	41.5"	45.5"	45.5"	36"
Commode opening	n/a	n/a	n/a	n/a	n/a	n/a	n/a	8"	8"	9"	n/a
Length	59.8"	62.8"	67.8"	72.3"	54.7"	60.5"	65.3"	54.7"	60.5"	65.3"	37"
Commode opening	n/a	n/a	n/a	n/a	n/a	n/a	n/a	11"	11"	13"	n/a
Back	29.2"	35.7"	41.7"	43.2"	n/a	n/a	n/a	n/a	n/a	n/a	n/a

Complete instructions for lift usage and sling procedures are found in the Invacare Operating Manual. Only Invacare slings and accessories should be used on Invacare lifts.

## Sling Binding Colors

Refer to chart below.



Note: Sling size and fit can vary significantly depending on patient weight and girth. These are general guidelines. Consult physician before sling selection.


Washing instructions are on each sling label. Washing temperature should not exceed 180°F.

Air-dry or dry at low temperatures.


All slings are constructed of polyester/-nylon fabrics which can be used in wet or dry environments.

## Sling Sizes

Color / Weight Chart

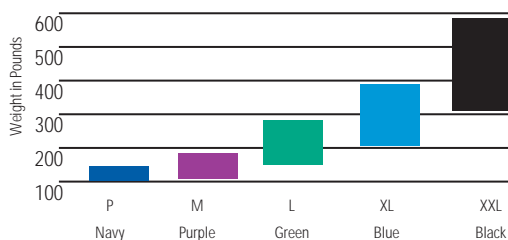

**Inspect sling before each use for wear, tears and loose stitching. Bleached, torn, cut, frayed or broken slings are unsafe and could result in injury. Discard immediately. Do not alter slings. Use only on Invacare lifts.**

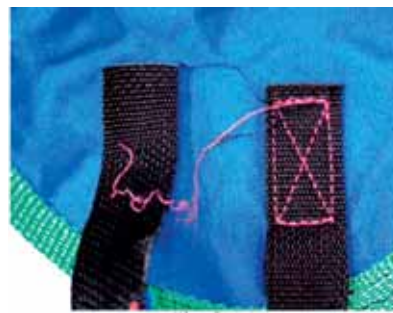

Use only Invacare lifts, slings and lift accessories to maintain patient safety. Use of other manufacturers products is prohibited and will void the warranty.

## Warning!

BEFORE using the Patient Lift, READ and UNDERSTAND the Owner's Manual for proper operation and safety procedures.

Use the sling that is recommended by the individual's doctor, nurse or medical attendant. Before lifting, check all sling straps for secure points of attachment on lift device. Do not exceed weight limitation posted on lift. Use only with genuine Invacare Patient Lifts.

Bleached, torn, cut, frayed or broken slings are unsafe and could result in injury. Discard immediately. Do not alter slings. Use only on Invacare lifts.

**For additional information on training manuals or videos, contact Invacare Corporation at (800)333-6900.**

**Invacare Corporation**  
www.invacare.com

### USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

### Canada

570 Matheson Blvd. E.  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2007 Invacare Corporation. All rights reserved.

Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted.

© 2007

Invacare Corporation  
Form No. 00-002  
Rev. 04/07



# Patient Sling Reference Guide



Yes, you can.®



- Full Body and Divided Leg Slings
- Toileting, Standing and Transport Slings
- Sizing Guide



Application: Use with Invacare® 9805 and Reliant 450 and 600 series power or manual lifts.  
 Models: 9805, RHA450-1, RHL450-1, RPA450-1, RPL450-1, RPL450-2, RPA600-1, RPL600-1, RPL600-2

## FULL BODY SLING

The Full Body Sling is an easy sling to use while transferring the resident from the bed to a wheelchair, toilet or shower-chair or from the floor to a bed. It offers total head and neck support and the padded lower section provides additional under-thigh comfort and support. Use the full body sling for residents who are totally dependent, partially-dependent, non-weight bearing, very heavy or have limited head control. Use this sling for amputee residents (in a reclined position). Choose from solid polyester fabric or bathing mesh, with or without commode opening. The mesh slings may be used for bathing or when the sling will be left under the resident. The solid polyester slings should not be left under the resident for extended periods because the padding retains body heat. Residents may be seated on "pressure-reduction" cushions. However, the efficacy of those cushions is diminished when solid padded slings are left under the resident.

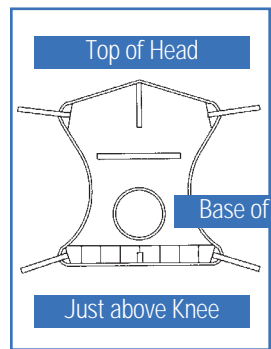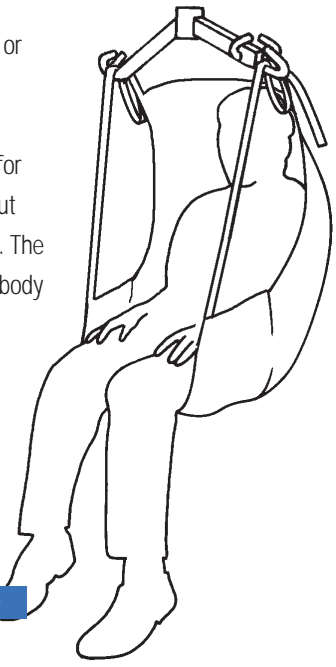


Model	Fabric & Style	Size/Trim Color
R110	Mesh, Full Body	Medium/Purple
R111	Mesh, Full Body	Large/Green
R112	Solid Fabric, Full Body	Medium/Purple
R113	Solid Fabric, Full Body	Large/Green
R114	Mesh, Full Body w/Opening	Medium/Purple
R115	Mesh, Full Body w/Opening	Large/Green
R116	Mesh w/Opening	Extra Large/Blue
R117	Solid Fabric	Extra Large/Blue
R140*	Heavy Duty Mesh	XXL/Black
R141*	Heavy Duty Mesh w/Opening	XXL/Black

\*Use Reliant 600 Lift


Model	Fabric & Style
9042	Solid Fabric, Full body
9043	Solid Fabric, Full Body w/Opening
9046	Mesh, Full Body
9047	Mesh, Full Body w/Opening

## DIVIDED LEG WITH HEADREST

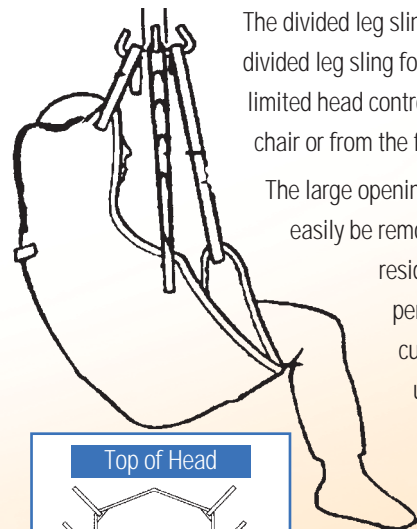


The divided leg sling is a secure, easy-to fit general purpose sling, designed for the majority of sling users. Use the divided leg sling for residents who are totally dependent, partially dependent, non-weight bearing, very heavy or have limited head control. Use this sling when transferring the resident from the bed to a wheelchair, recliner or shower-chair or from the floor to a bed.

The large opening also accommodates toileting. This versatile sling provides full head and neck support and can easily be removed or applied to a resident in a seated or supine position. The solid fabric has extra padding for residents with delicate skin. The solid polyester slings should not be left under the resident for extended periods because the padding retains body heat. Residents may be seated on "pressure-reduction" cushions. However, the efficacy of those cushions is diminished when solid padded slings are left under the resident.

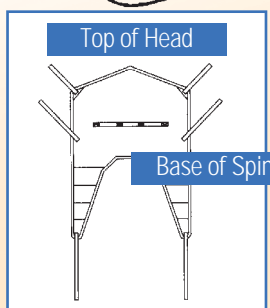


Model	Fabric & Style	Size/Trim Color
R100	Solid Fabric, Divided Leg w/Headrest	Medium/Purple
R101	Solid Fabric, Divided Leg w/Headrest	Large/Green
R102	Solid Fabric, Divided Leg w/Headrest	Extra Large/Blue
R100P	Solid Fabric, Divided Leg w/Headrest	Petite/Navy

## TOILETING SLING

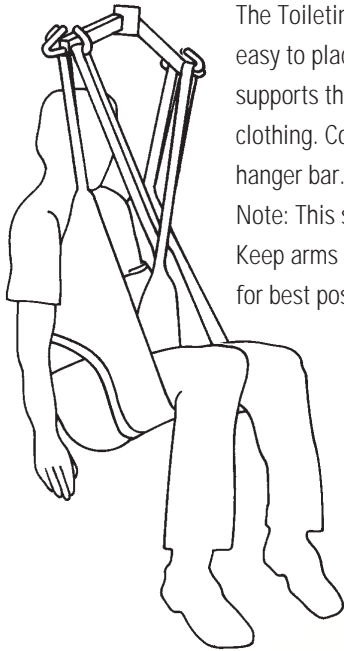

The Toiletting Sling offers superb access for quick toileting and hygiene. It is easy to place and remove in a bed or chair situation. The fully padded design supports the back, offers under thigh comfort and allows for removal of clothing. Color coded positioning straps help align the sling properly on the hanger bar.

Note: This sling requires the patient to have FULL head and neck control. Keep arms on outside of toiletting sling. Place leg sections 3-4" above knee for best positioning.

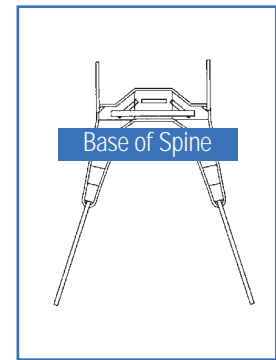


Model	Fabric & Style	Size/Trim Color
R121	Solid Fabric, Toileting	Large/Green

**Application:** Use only with Reliant Stand Up Lift.

**Models:** RPS350-1, RPS350-2, RPS440-1

## STANDING SLING R130

The R130 and R131 standing slings are for residents who are partially dependent, have at least 60% weight bearing capacity, have head and neck control, are cooperative, can sit up on the edge of the bed (with or without assistance), and are able to bend at the hip, knees and ankles. Use the R130 or R131 standing slings for standing assistance, quick toileting, weight bearing practice and for transferring resident from bed to chair, or chair to bed.

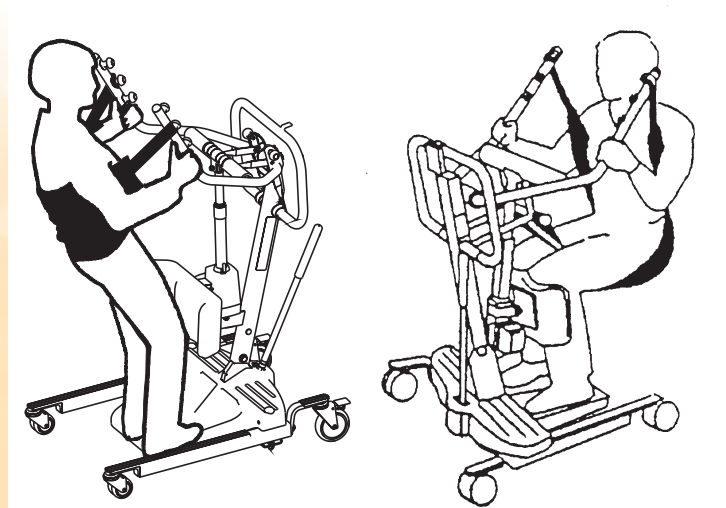


## TRANSPORT SLING R134

The R134 and R136 transport slings are for transferring residents with minimal weight-bearing capacity, who have head and neck control, are cooperative, can sit up on the edge of the bed (with or without assistance), and are able to bend at the hips, knees and ankles. Use the R134 and R136 transport sling for transferring resident from bed to chair, chair to car or general in-room transport. Do not raise resident to a full standing position while using the transport sling.

## SLING POSITIONING (Reclined)

Place fan-folded or rolled sling under patient with grab strap away from skin. Roll patient from side to side, smoothing each half of the sling. Bed rail on opposite side should be raised. Refer to Invacare Sling Owner's Manual for detailed instructions.

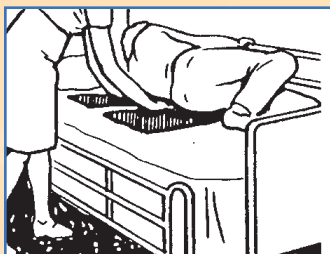

Full Body (above)



Divided Leg (right)

## SLING POSITIONING (Seated)

For Divided Leg or Toiletting Sling, slide behind seated patient until belt is at waist level or bottom of back section is at the base of the spine. Place leg sections under thighs, smoothing fabric. Attach to hanger bar in one of the three methods: crossed between legs, under both legs or single loop under each leg.

